

WORK PERMIT

Department of Chemical and Biological Engineering

化學及生物工程學系

Project Title : Persistence Study of Staphylococcus aureus
towards Antibiotics by Shotgun Proteomics

Researcher(s) : Jordy Evan SULAIMAN

Supervisor(s) : Prof. Henry H. N. Lam

Work Plan No. : 17057

Date of Approval : Aug 15th 2017

Date of Revalidation : N/A

Signature of Approval : *Pauline Leung*

Mrs Pauline S T LEUNG
DSO



NEW

The Hong Kong University of Science & Technology
Department of Chemical and Biomolecular Engineering

Persistence Study of
***Staphylococcus aureus* towards**
Antibiotics by Shotgun Proteomics

Work Plan #17057

Researcher: Jordy Evan SULAIMAN

Supervisor: Prof. Henry H. N. Lam

July, 2017

1. General Information

Name of Researcher:	SULAIMAN, Jordy Evan
Name of Project Supervisors:	Prof. Henry H. N. Lam
Project Title:	Persistence Study of <i>Staphylococcus aureus</i> towards Antibiotics by Shotgun Proteomics
Research Area:	Proteomics
Proposed Start Date:	August 18, 2017
Location:	Room 7110 and Bioengineering Laboratory

2. Experiment/Project Description

Methicillin-resistant *Staphylococcus aureus* (MRSA) has become a global health threat causing fatal bacterial infections in hospitals. However, until today, the response mechanism of *S. aureus* towards antibiotics, especially β -lactams antibiotics, remains unclear. In order to treat MRSA infection, the antibiotic resistance mechanism and defense in *S. aureus* need to be studied further. In this project, quantitative proteomics will be employed to study the antibiotic responses in *S. aureus*. Furthermore, investigation regarding whether the persister cells survive antibiotic treatment by reducing metabolism or due to active response to the antibiotic stress will be done.

In this work, two strains of *S. aureus* (ATCC43300/ MRSA and ATCC25923/ MSSA) will be cultured, and the cell pellets will be obtained. The cell pellets will be lysed to extract the proteins, and the proteomes will be analyzed by standard shotgun proteomics techniques (LC-MS/MS). The method has been used in the research field for investigating the proteome changes of cells, and the knowledge obtained from the experiment can provide new insights into the major factors in the antibiotic response of *S. aureus*.

After the mass spectrum is successfully obtained, it will be compared with the protein spectrum existed in the Spectra Library. Comparing experimental mass spectra with Spectra Library require computational calculation, and the protein produced by *S*

.aureus in response to antibiotics can be identified. Furthermore, the proteins will also be quantified by spectral counting.

3. Equipment List

1. Biochemical Incubator	2. Incubator shaker
3. Biological Safety Cabinet	4. Chemical Fume Hood
5. Centrifuges	6. SpeedVac
7. Micro-plate Reader	8. Ultrasonic Homogenizer
9. High Performance Liquid Chromatography (HPLC)	10. Linear Trap Quadrupole (LTQ) Mass Spectrometry
11. Quadrupole Time-of-Flight (Q-TOF) Mass Spectrometry	

4. Experimental Procedures

4.1 Cell Culture

4.1.1 Autoclave all glassware that will be used for experiments under 121°C for 20 minutes (103.4 kPa).

4.1.2 Preparation of Mueller-Hinton liquid media

Reagent	Amount to add (for 100 mL)
Acid hydrolysate of casein	1.75 g
Beef extract	0.3 g
Starch	0.15 g
H ₂ O	100 mL

4.1.3 *Staphylococcus aureus* cell culture

Pick a single *S. aureus* colony from a normal agar plate and inoculate a starter culture of 2-5 mL Mueller-Hinton liquid media with different antibiotics in different test-tubes. Incubate for approximately 8-12 h at 37 °C with vigorous shaking in nearly 300rpm.

4.2 Protein Extractions from bacteria

Experiment performed inside BSC II

4.2.1 Harvest the *S. aureus* cells by centrifugation at $6000 \times g$ for 15 minutes under 4°C and aspirate the supernatant. Suctioning the supernatant must be conducted inside the BSC II to make sure no *S.aureus* contamination to outside environment. For *S. aureus* preparation inside the BSC, always wear mask to prevent inhalation of *S. aureus*.

4.2.2 Resuspend the pellet in cold PBS buffer and centrifuge at $6000 \times g$ for 10 minutes under 4°C and aspirate the supernatant. Repeat 2 times.

4.2.3 Add ice-cold cell lysis buffer (8 M urea) and resuspend the pellet. Incubate on ice for 10 minutes. Close the cap and make sure there is no spillage on outside of the micro centrifuge tubes. The micro centrifuge tube must be sterile before taken out from the BSC.

Experiment performed on the benchtop outside BSC

4.2.4. Put the *S. aureus* pellet sample (in the closed micro centrifuge tube) inside liquid N₂. Repeat 2 times.

4.2.5 Vortex tubes briefly and proceeds to sonication for 10 min at 4°C.

4.2.6 Transfer the supernatant to a new tube and centrifuge samples at $6000 \times g$ for 15 minutes at 4°C to remove any insoluble material.

4.2.7 Add ice-cold acetone to precipitate the proteins.

4.2.8 Aspirate the supernatant and redissolve the proteins in buffer (4 M urea and 30 mM Tris-HCl, pH 6.5).

4.2.9 Take an aliquot for the protein quantification.

4.2.10 Repeat step 4.2.1 to 4.2.8 for other antibiotic samples.

4.3 Sample preparation

4.3.1 Reduce the proteins with dithiothreitol (DTT, 10 mM final concentration) at 37°C for 3 hours.

4.3.2 Alkylate the proteins with iodoacetamide (IAA, 20 mM final concentration) in dark for 1 hr. The alkylation reaction is quenched by adding DTT (10 mM final concentration) again.

4.3.3 Dilute the samples to a concentration of urea less than 1 M.

4.3.4 Digest the proteins with sequencing grade modified trypsin (1: 50, w/w) at 37°C overnight.

4.3.5 Acidify samples with 10% formic acid to a final concentration of 0.5% (v/v).

4.3.6 Desalt the sample with C18 reverse-phase ZipTip.

4.3.7 Dry the sample with SpeedVac and store it at -20°C before use.

4.4 Analysis by mass spectrometry (LC-MS/MS)

4.4.1 Run the samples into Thermo scientific LTQ Velos platform mass spectrometer which is interfaced to a nanoelectrospray ion source coupled to a Thermo Accela LC.

4.5 Autoclaving waste containing *S.aureus*

4.5.1 Autoclave waste containing *S.aureus* bacteria under 121°C for 3 hours (103.4 kPa).

5. Procedure Template

Experimental Procedure No.	Experimental Procedure Description	Scale (Mass/Volume)	Location (Fumehood, benchtop, etc.)	Method (New or Existing)
4.1.2	Preparation of Mueller-Hinton liquid media	Acid hydrolysate of casein = 1.75 g Beef extract = 0.3 g Starch = 0.15 g H ₂ O = 100 mL	Room 6104, benchtop	Existing
4.1.3	<i>S. aureus</i> cell culturing	2~5 mL	Room 7110, BSC	Existing
4.2.1	Harvesting <i>S. aureus</i>	-	Room 7110, BSC	Existing
4.2.2 4.2.3	Lysing and Killing <i>S. aureus</i> cells	PBS buffer-5 mL Lysis buffer-2 mL	Room 7110, BSC	Existing
4.2.7	Protein extraction	Acetone-5 mL	Room 7110, benchtop	Existing
4.2.9	Protein quantification	800 µL	Room 6104, Microplate reader	Existing
4.3	Sample preparation	DTT-2 µL IAA-2 µL Sequencing grade modified trypsin-2 µL Formic acid-2 µL	Room 7110, benchtop	Existing
4.4	Analysis by mass spectrometry (LC-MS/MS)	N/A	Room 7101 and BioCRF	Existing

6. Hazard and Operability Analysis (HAZOP)

Hazards and Operability Analysis							
4.1 Cell Culture							
No.	Hazards	Hazards Effect	Severity	Probability	Risk	Minimize Risk By	Residual Risk
1	Contact with chemicals	Causes sever skin and eye burns	H	M	H	Wear protective gloves, face shield and lab coats; Conduct experiments in the fume hood	L
2	Autoclave	High pressure steam burns	H	M	H	Strictly follow the safety requests and do not overload the autoclave	L
3	Contact with bacteria	Causes illness and infection	H	M	H	Wear protective gloves, face shield and lab coats; Conduct experiments in the biosafety cabinet	L
FINAL ASSESSMENT:						OVERALL RISK:	L
4.2 Protein Extraction							
1	Contact with chemicals	Causes sever skin and eye burns	H	M	H	Wear protective gloves, face shield and lab coats; Conduct experiments in the fume hood	L
2	Centrifuge	Rotor Failure: the rotor that breaks can spin out of control and hit laboratory personnel.	H	M	H	Strictly follow the safety requests.	L

		Spillage: the tubes may spray over the machine parts.				Balance all samples as closely as possible Set the rotors under the maximum speed. OVERALL RISK: L
FINAL ASSESSMENT:						
4.3 Sample Preparation						
1	Contact with chemicals	Causes sever skin and eye burns	H	M	H	Wear protective gloves, face shield and lab coats; Conduct experiments in the fume hood L
2	SpeedVac	Rotor Failure: the rotor that breaks can spin out of control and hit laboratory personnel. Spillage: the tubes may spray over the machine parts.	H	M	H	Strictly follow the safety requests. Balance all samples as closely as possible Set the rotors under the maximum speed. L
3	Autoclave	High pressure steam burns	H	M	H	Strictly follow the safety requests and do not overload the autoclave L
FINAL ASSESSMENT:						
OVERALL RISK: L						

Severity: L=Low (Minor injuries, first aid), M=Middle (Hospitalization, medical leave), H=High (Serious injuries, fatality).

Probability: L=Low (Unlikely), M=Middle (Possible), H=High (Very likely).

Risk=Severity × Probability, the product follows the higher severity or probability

7. Operating Conditions

- Autoclaving equipment: 121°C, 103.4 kPa, 20 min
- Autoclaving waste: 121°C, 103.4 kPa, 3 hours
- Cell Culture: 37°C, shaking, 8~12 hr
- Pellet Cells: 6000 × g, 4°C, 10~15min
- Reduce Proteins: Dithiothreitol (DTT) (10 mM final concentration), 37°C, 3 hours
- Alkylate Proteins: Iodoacetamide (IAA) (20 mM final concentration), in dark, 1 hr
- Digest Proteins: 37°C, overnight

8. Services List

- Electricity (AC 220V, 50Hz)
- Tap water
- Double Deionized Water
- Chemical Fume Hood
- Biological Safety Cabinet (Level 2)
- Linear Trap Quadrupole (LTQ) Mass Spectrometry

9. Chemicals List

Chemical	Purity	Quantity per Experiment
Acid hydrolysate of casein	99.999%	1.75 g
Beef extract	99%	0.3 g
Starch	99.99%	0.15 g
Urea	95%	50 µL (4M solution)
Tris-HCl	99%	10 mL (30 mM solution)
Dithiothreitol	99.5%	1 mg
Iodoacetamide	99%	1 mg
Formic Acid (CH ₂ O ₂)	98%	1 mL
Acetone (C ₃ H ₆ O)	99.9%	5 mL

Trypsin	2,500 USP units/mg	3 mg
Phosphate Buffered Saline (PBS)	99.9%	10 mL
Oxacillin Sodium Salt (C ₁₉ H ₁₈ N ₃ NaO ₅ S)	99%	200 µg
Ampicillin sodium salt (C ₁₆ H ₁₈ N ₃ NaO ₄ S)	99%	200 µg
Amoxicillin trihydrate (C ₁₆ H ₁₉ N ₃ O ₅ S · 3H ₂ O)	99%	200 µg
Cephalexin Hydrate (C ₁₆ H ₁₇ N ₃ O ₄ S · xH ₂ O)	99%	200 µg
Sulfathiazole Sodium salt (C ₉ H ₈ N ₃ NaO ₂ S ₂)	99%	200 µg
Tetracycline hydrochloride (C ₂₂ H ₂₄ N ₂ O ₈ · HCl)	99%	200 µg
Erythromycin (C ₃₇ H ₆₇ NO ₁₃)	99%	200 µg
Levofloxacin (C ₁₈ H ₂₀ FN ₃ O ₄)	99%	200 µg

10. BIOLOGICAL AGENTS

BIOLOGICAL AGENT	BIOLOGICAL SAFETY LEVEL	REFERENCE FOR BSL LEVEL (Attachment: Compulsory)
<i>Staphylococcus aureus</i> ATCC25923 (MSSA)	2 (Pathogenic, not airborne)	https://www.atcc.org/Products/All/25923.aspx
<i>Staphylococcus aureus</i> ATCC43300 (MRSA)	2 (Pathogenic, not airborne)	https://www.alliance-bio-expertise.com/eshop/staphylococcus-aureus-subsp-aureus-atcc-43300-kwik-stik-x2.html

See attached SDS for *S. aureus* for a complete assessment.

(NOTE: If reputable bodies differ in their BSL Level for any Biological Agent, the Highest BSL must be used for evaluating risks and deciding on which safety procedures are to be enforced.)

11. Summary of Relevant Hazards and Incompatibilities

Chemical	Summary of Hazards	Incompatibilities
Acid hydrolysate of casein	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Beef extract	This substance is not classified as dangerous according to Directive 67/548/EEC.	Strong oxidizing agents
Starch	This substance is not classified as dangerous according to Directive 67/548/EEC.	Strong oxidizing agents
Urea	This substance is not classified as dangerous according to Directive 67/548/EEC.	Strong oxidizing agents
Tris-HCl	This substance is not classified as dangerous according to Directive 67/548/EEC.	Bases, Oxidizing agents
Dithiothreitol	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Bases, Oxidizing agents, Reducing agents, Alkali metals
Iodoacetamide	Toxic if swallowed. May cause sensitization by inhalation and skin contact.	Strong acids, Strong bases, Strong oxidizing agents, Strong reducing agents
Formic Acid (CH₂O₂)	Flammable liquid and vapor. Causes severe skin burns and eye damage.	Strong oxidizing agents, Strong bases, Powered metals
Acetone (C₃H₆O)	Highly flammable. Irritating to eyes. Repeated exposure may cause skin dryness or cracking. Vapors may cause drowsiness and dizziness.	Bases, Oxidizing agents, Reducing agents, Acetone reacts violently with phosphorous oxychloride
Trypsin	May cause allergy or asthma symptoms or breathing difficulties if inhaled. Irritating to skin and eyes	N/A
Phosphate Buffered Saline (PBS)	Not a hazardous substance or mixture according to EC-	Strong oxidizing agents, Strong acids

	directives 67/548/EEC or 1999/45/EC.	
Oxacillin Sodium Salt (C ₁₉ H ₁₈ N ₃ NaO ₅ S)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Ampicillin sodium salt (C ₁₆ H ₁₈ N ₃ NaO ₄ S)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Oxidizing agents
Amoxicillin trihydrate (C ₁₆ H ₁₉ N ₃ O ₅ S · 3H ₂ O)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Cephalexin Hydrate (C ₁₆ H ₁₇ N ₃ O ₄ S · xH ₂ O)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Sulfathiazole Sodium salt (C ₉ H ₈ N ₃ NaO ₂ S ₂)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Tetracycline hydrochloride (C ₂₂ H ₂₄ N ₂ O ₈ · HCl)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents
Erythromycin (C ₃₇ H ₆₇ NO ₁₃)	This substance is not classified as dangerous according to Directive 67/548/EEC.	Strong oxidizing agents
Levofloxacin (C ₁₈ H ₂₀ FN ₃ O ₄)	Harmful if swallowed. Irritating to eyes, respiratory system and skin.	Strong oxidizing agents

12. Waste List

1. Cell culture (liquid) containing Mueller-Hinton liquid medium, autoclaved and dispose to the sink.
2. Used PBS (liquid), autoclaved and disposed in the sink.
3. Used acetone (liquid), disposed in the organic waste container of non-halogenated solvent.
4. Used pipette and microcentrifuge tubes (solid), disposed in the waste container of biological waste and disposed as domestic waste after autoclave

13. Assessment of Significant Risks

May exposed to toxic or irritating chemicals. Goggles and gloves should be always put on.

14. Safety Precautions

- Personal Protective Equipment
 - A lab coat, safety glasses and respirator must be worn throughout the experiment. Rubber glove should be worn according to the nature of chemicals being handled.
- Fume Extraction System is required and in Place
 - The reaction system is run in the fume cupboard throughout the course of the experiment.
- Warning sign required
 - Chemicals have to be clearly labeled, e.g. harmful, toxic, flammable, and irritant.
- Safety Training Required
 - Safety training courses, included Chemical Safety (I), (II) and Biological Safety organized by the HSEO should be attended.
- Specific training required
 - None
- Emergency shutdown procedures
 - Follow the emergency shutdown procedures, absorb any liquid spillage using the laboratory caustic spill kit, transfer used materials to an appropriately labeled container and inform the person designated responsible for laboratory.
- If connected with chemicals and reagents
 - Immediately flush eyes or skin with top water at least 15 minutes and take off contaminated clothing as far as possible.

15. Action in Case of Abnormal or Emergency Situations

15.1 Service Failure

Close the main valve on the air cylinder and reactant inlet.

15.2 Action in case of fire or explosion

(1) When hear the fire alarm

- Remain calm and check if there is any sign of fire in the vicinity.
- If you see fire or smoke, or hear the announcement asking you to evacuate, follow the evacuation procedures below.
- If there is no sign of a fire, stay alert and pay attention to announcement until the fire alarm is silenced.
- Evacuate if the alarm has sounded for more than two minutes.

- If you hear the buzzer sound which indicates fire alarm is activated in an adjacent fire zone, stay alert and pay attention to announcement.

- If you hear both the buzzer and the fire alarm, you are near the boundary of fire zones, treat as if you hear the fire alarm.

(2) If discover a fire

- Activate the fire alarm by pressing the break glass fire alarm button.
- Report to Security Control Centre by dialing 8999.
- Alert other people. If SAFE to do so, try to put out the fire by firefighting equipment.

- DO NOT take any personal risk. If the fire gets beyond your control, evacuate immediately by following the procedures below.

- Close the door of the room on fire.

(3) When there is a fire and need to evacuate

- Remain calm. Walk, do not run, especially when travelling on staircases.
- Immediately leave the building and go to the assembly point using the nearest exit.

- Try to help those who may have difficulties travelling such as disabled and pregnant persons.

- DO NOT USE THE LIFTS.

- Report to your Fire & Safety Officer at the assembly point as far as practicable.

- Do not return to the building until permission is given by the Fire Services Department Officer in charge of the scene.

(4) Fire Fighting Equipment

- Water from the hose reels is good for wood, paper and structural fire, but NOT for oil, electrical or metal fire.

- The most common fire extinguishers on campus are the carbon dioxide types (black containers) which are good for general purposes, including oil and/or electrical fire.

- Some laboratories have dry powder fire extinguishers (blue container), which are good for chemical and/or metal fire.

- Sand (lab buckets) can be used to contain flammable liquid as well as put out a fire, including metal fire.

- Fire blanket can be used when someone's clothing catches fire.

15.3 Action in case of hazardous chemical spill in a research laboratory

- Alert co-workers.

- If safe to do so, (1). Confine the spill with appropriate materials. (2). Turn off remotely all heat/ignition sources if flammable vapor is involved.

- Ask for assistance if necessary.

- Press the Emergency Ventilation button (do not activate this button in case of fire).

- Inform the Security Control Centre by dialing 8999 from a safe location.

- Evacuate everyone in the affected area. Leave contaminated clothing and close the door.
- Activate local warning system to prevent others from entering the room.
- If possible, maintain a safe distance from the scene, keep the entrance or access routes in sight and help to prevent entry to the affected room.
- If conditions allow, remain to assist the emergency response team.

15.4 Action in case of other abnormal situations

In case of accidental contact with chemicals, flush eyes or skin with copious amount of water immediately for at least 15 minutes while removing contaminated clothing and shoes. For inhalation, remove to fresh air. If swallowed, wash mouth with water. Contact nearest physician.

CBME Risk Assessment Audit Declaration

I, SULAIMAN Jordy Evan, have read and understood the relevant sections of the HKUST Safety and Environmental Protection manual. I am aware of my responsibilities with regard to the health and safety of myself and others and state that, to the best of my knowledge, the information provided in this work plan is complete and correct:

I have attended and passed the following safety courses :

- ☒ 1. Hazardous Waste Management (**mandatory**) __MC03
- ☒ 2. Chemical Safety for Laboratory Users (**mandatory**) __MC07
- ☒ 3. Fire Safety (**mandatory**)
- ☐ 4. Laser Safety __MC04
- ☒ 5. Biological Safety __MC06
- ☐ 6. Pressure Safety __MC05
- ☐ 7. Others : _____

Are any of the following categories relevant to the proposed research?

- | | | |
|--|---|--|
| <input type="checkbox"/> High voltage power supplies | <input type="checkbox"/> Yes | <input checked="" type="checkbox"/> No |
| <input type="checkbox"/> Biologically active materials | <input checked="" type="checkbox"/> Yes | <input type="checkbox"/> No |
| <input type="checkbox"/> Radioactive materials or ionising radiation sources | <input type="checkbox"/> Yes | <input checked="" type="checkbox"/> No |
| <input type="checkbox"/> Non-ionising radiation (UV, microwaves, lasers) | <input checked="" type="checkbox"/> Yes | <input type="checkbox"/> No |
| <input type="checkbox"/> Highly toxic, carcinogenic or mutagenic materials | <input type="checkbox"/> Yes | <input checked="" type="checkbox"/> No |
| <input type="checkbox"/> Highly flammable or explosive materials | <input type="checkbox"/> Yes | <input checked="" type="checkbox"/> No |
| <input type="checkbox"/> Operation under extreme pressure or temperature | <input type="checkbox"/> Yes | <input checked="" type="checkbox"/> No |


Status of Researcher

1. ☒ HKUST **CBME** ☐ HKUST _____ ☐ Non HKUST _____ ☐ Others _____
2. ☐ Ph.D. ☒ M.Phil. ☐ M.Sc.
☐ FYP ☐ UROP ☐ Student Helper
☐ Visiting Scholar ☐ Research Assistant ☐ Research Associate
☐ Intern-PG ☐ Exchange-PG ☐ Post Doc
☐ Intern-UG ☐ Exchange-UG ☐ Faculty
☐ UG Research
☐ _____

Insurance Coverage

All researchers *other than undergraduate or postgraduate students* must provide a copy of their letter of appointment to ensure that they have official status within the Department and that they have appropriate insurance coverage. Any researchers employed by anyone other than the Department of Chemical and Biomolecular Engineering at HKUST, must provide full details of their insurance coverage by their employer.

Assessment checked by Technical staff associated with ALL laboratories in which the planned EXPERIMENTAL work will be undertaken.

Laboratory	Responsible Technician (Name)	Signature	Date
7110	Kam Tim TANG		

Please list ALL laboratories in which ANY ANALYTICAL work will be undertaken:

Bioengineering Laboratory (Room 6104),


BioCRF (Room 6127)

Laboratory	Responsible Technician Name	Equipment	Training Required Y/N	Training Completed Y/N
6104	Inez S M Tsui	BCA assay	Y	N
6127	Joyce Wong	LC-MS	Y	N


* **Note:** After completion of training, a record of training completion will be recorded and kept by the relevant technician

I have had my work plan checked by the relevant technical staff associated with all laboratories in which the planned work will be undertaken and also by my Faculty research supervisor. I have considered their suggestions and feedback and have modified the work plan accordingly.

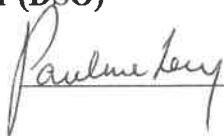
Researcher

Name	Signature	Date
<u>SULAIMAN, Jordy Evan</u>	<u></u>	<u>26 July 2017</u>

Assessment checked by research supervisor(s)

Name	Signature	Date
<u>Prof. Henry H. N. Lam</u>	<u></u>	<u>15 Aug 2017</u>

Authorized by CBME Departmental Safety Officer (DSO)

Name	Signature	Date
<u>Mrs. Pauline S.T. Leung</u>	<u></u>	<u>Aug 15th 2017</u>

PATHOGEN SAFETY DATA SHEET

Staphylococcus aureus

CHARACTERISTICS

	Gram-positive cocci, usually occurs in clusters, non-spore forming, non-motile, coagulase positive, facultative anaerobes.
<i>Morphology</i>	
<i>Growth Conditions</i>	Tryptic Soy Broth

HEALTH HAZARDS

<i>Host Range</i>	Humans and Animals.
<i>Modes of Transmission</i>	Ingestion of food containing enterotoxins, contact with nasal carriers, contact with draining lesions or purulent discharges, also spread by person-to-person contact; Indirectly by contact with fomites, indirectly or directly by contact with infected animals.
<i>Signs and Symptoms</i>	<i>Accidental ingestion:</i> Violent onset of severe nausea, cramps, vomiting, and diarrhea if preformed enterotoxin is present. <i>Surface infections:</i> Impetigo, folliculitis, abscesses, boils, infected lacerations. <i>Systemic infections:</i> onset of fever, headache, myalgia, can progress to endocarditis, meningitis, septic arthritis, pneumonia, osteomyelitis, sepsis.
<i>Infectious Dose</i>	Virulence varies for different strains.
<i>Incubation Period</i>	30 minutes to 8 hours when consuming contaminated food with enterotoxin. Otherwise, typically 4 to 10 days. Disease may not occur until several months after colonization of mucosal surfaces.

MEDICAL PRECAUTIONS/TREATMENT

<i>Prophylaxis</i>	Hand-hygiene; Elimination of nasal carriage by using topical mupirocin. Mupirocin also eliminates transient hand carriage by eliminating the mucosal reservoir.
<i>Vaccines</i>	None
<i>Treatment</i>	Incision and drainage for localized skin infections; antibiotic therapy for severe infections; Many strains resistant to antibiotics; Sensitivity must be determined for each strain.
<i>Surveillance</i>	Monitor for signs of food poisoning when ingestion occurs. Monitor for skin inflammation; isolation of organism from wound, blood, CSF or urine.
<i>MSU Requirements</i>	Report any exposures

LABORATORY HAZARDS

<i>Laboratory Acquired Infections (LAIs)</i>	29 reported cases up to 1973 with 1 death. Most common cause of laboratory infection was accidental self-exposure via the mucous membranes by touching contaminated hands to face or eyes.
<i>Sources</i>	Contaminated food, blood, abscesses, lesion exudates, CFS, respiratory specimen, feces, and urine

SUPPLEMENTAL REFERENCES

<i>Canadian MSDS:</i>	http://www.phac-aspc.gc.ca/lab-bio/res/psds-ftss/index-eng.php
<i>BMBL:5th Edition</i>	http://www.cdc.gov/biosafety/publications/bmbl5/BM5.pdf
<i>CDC</i>	http://www.cdc.gov/HAI/organisms/staph.html

CONTAINMENT REQUIREMENTS

<i>BSL2</i>	For all procedures involving known or potentially infected cultures.
<i>ABSL2</i>	For all procedures utilizing infected animals

SPILL PROCEDURES

<i>Small</i>	Notify others working in the lab. Remove and don new PPE. Cover area of the spill with absorbent material and add 10 % Bleach. Allow 30 minutes hour of contact time. After 30 minutes and then cleanup and dispose of materials.
<i>Large</i>	For assistance, contact MSU's Biosafety Officer (406-994-6998) or Safety and Risk Management (406-994-2711).

EXPOSURE PROCEDURES

<i>Mucous membrane</i>	Flush eyes, mouth or nose for 5 minutes at eyewash station.
<i>Other Exposures</i>	Wash area with soap and water for 5 minutes.
<i>Reporting</i>	Immediately report incident to supervisor, complete a first report of injury report, and submit to Safety and Risk Management. During business hours: Montana Occupational Health 2075 Charlotte St. Suite 3 Bozeman, MT After business hours: Bozeman Deaconess Hospital Emergency Room 915 Highland Blvd Bozeman, MT
<i>Medical Follow-up</i>	

VIABILITY

<i>Disinfection</i>	Susceptible to 10 % Bleach, 70 % ethanol, and 2 % glutaraldehyde, chlorohexadine, formaldehyde, and 0.25 % benzalkonium chloride.
<i>Inactivation</i>	Inactivated by dry heat (1 hour at 160-170°C).
<i>Survival Outside Host</i>	Carcass and organs – 42 days; Skin – 30 minutes to 38 days; meat products – 60 days; floor – less than 7 days; glassware – 46 hours; sunlight – 17 days; UV light – 7 hours.

PERSONAL PROTECTIVE EQUIPMENT (PPE)

<i>Minimum PPE Requirements</i>	At minimum, gloves, closed toed shoes, lab coat, and appropriate face and eye protection prior to working with <i>S. aureus</i> . Additional PPE may be required depending on lab specific SOPs.
<i>Additional Precautions</i>	Avoid injuries from contaminated sharp instruments, bites and scratches from infected animals, and direct contact with open skin or lesions of skin.

